


# Organophosphorus pesticides

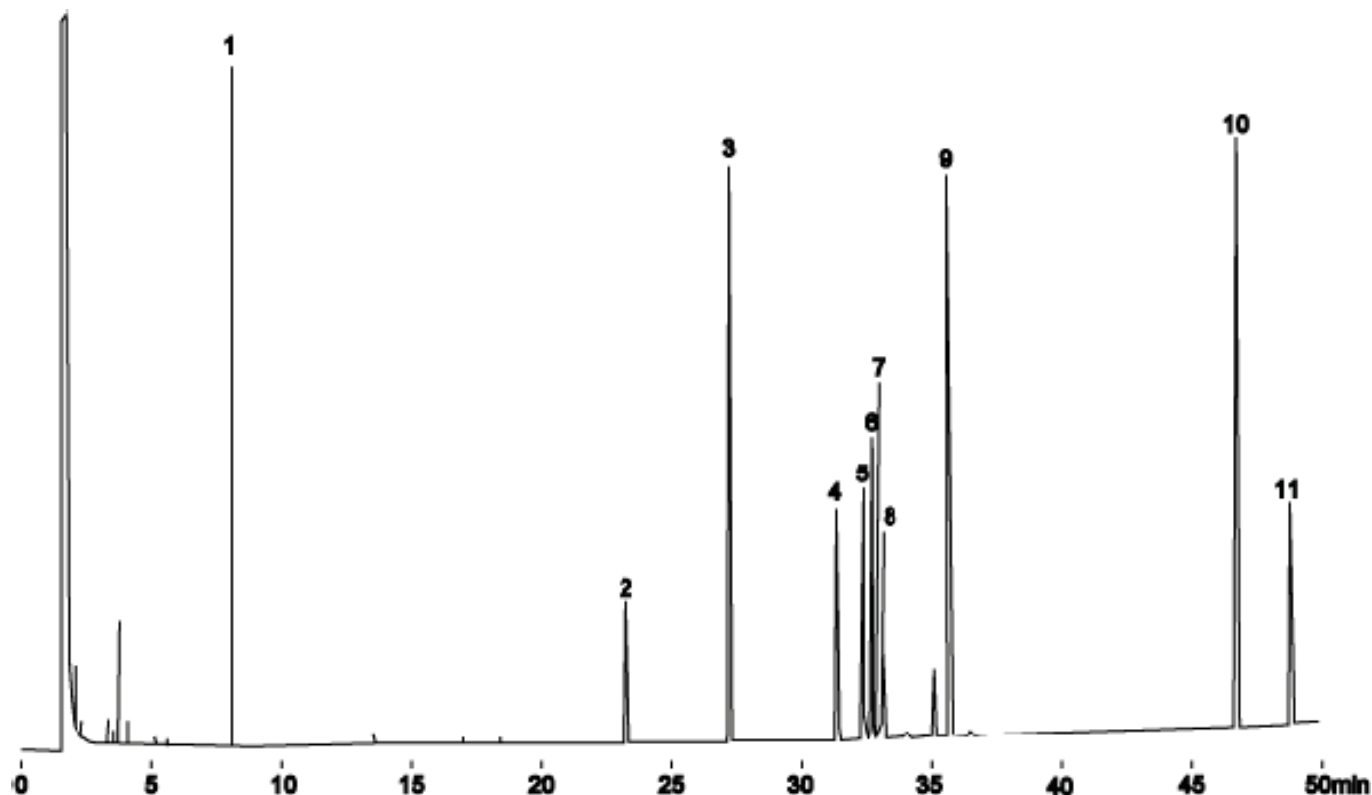

System : GC/FID  
 Column : InertCap 1  
 0.25mmI.D. × 30m  
 df=0.4μm  
 Column Temp. : 100°C→3°C/min→270°C  
 Carrier Gas : He 80kPa  
 Injection : Split 1 : 50 1μl  
 Detection : FID Range 1 × 256  
 Sample : Pesticides

Solutes : 1. Dichlorovos  
 2. Dimethoate  
 3. Diazinon  
 4. Fenitrothion  
 5. Malathion  
 6. Fenthion  
 7. Parathion  
 8. Chlorpyriphos  
 9. Chlorfenvinphos (CVP)  
 10. Phenthoate (PAP)  
 11. EPN  
 12. Phosalone