


Flavor by Japanese Standards of Food Additives

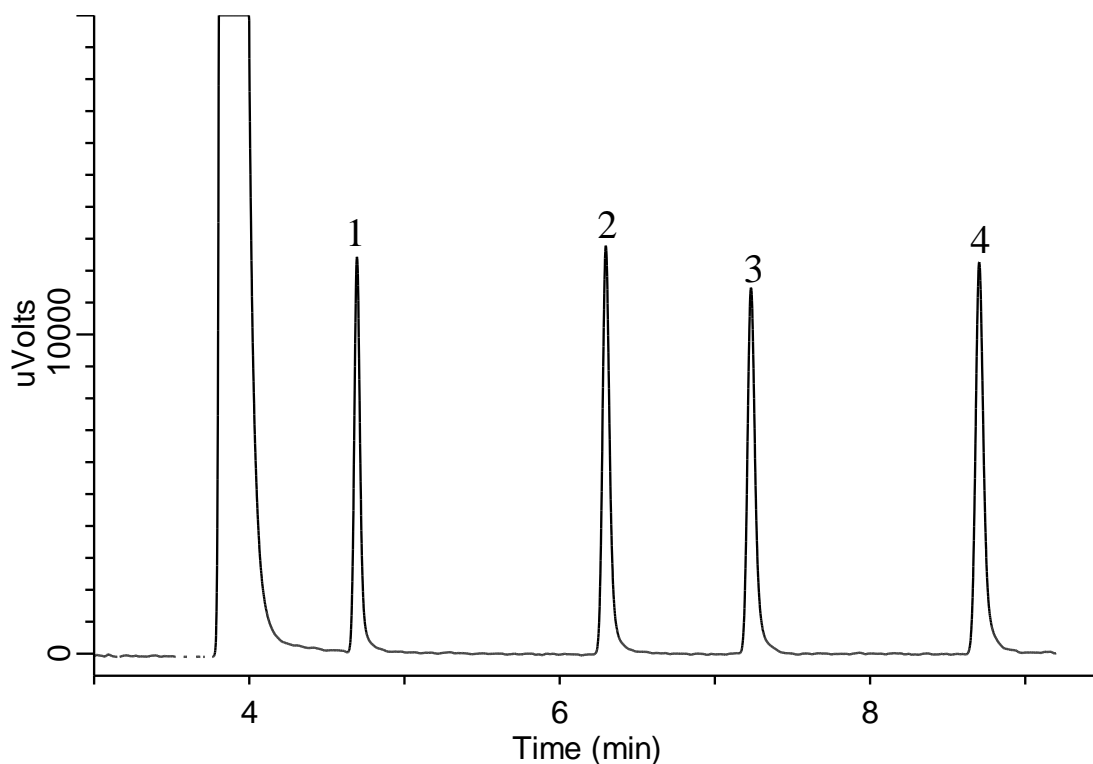

Conditions

System : GC-FID
Column : InertCap Pure-WAX
 0.25 mm I.D. × 30 m df = 0.25 μm
Col. Cat. No. : 1010-68142
Col. Temp. : 50 °C(5 min hold) - 5 °C/min - 230 °C
Carrier Gas : He 70 kPa
Injection : Split flow 96 mL/min
 250 °C
Detection : FID Auto Range
 250 °C
Sample Size : 1000 μg/mL in Methanol
 1 μL

Analyte

1. Ethyl propionate
2. Ethyl butyrate
3. Butyl acetate
4. Isoamyl acetate