


Vitamin E (Tocopherols) in Margarine

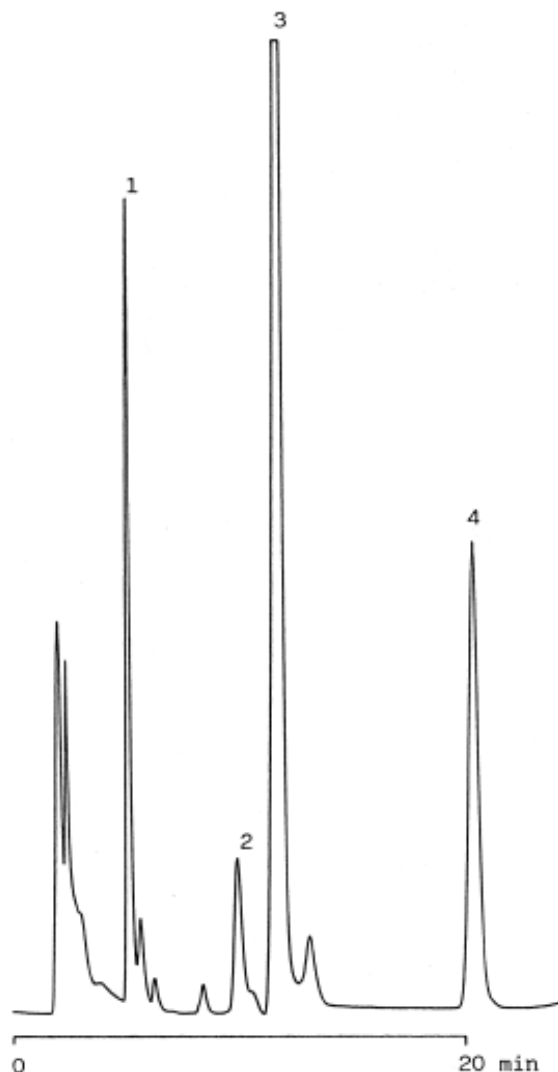

Column : Inertsil SIL-150A
5µm
5020-01011

Column Cat. No.:

Eluent : n-Hexane/Ethanol = 99.7/0.3

Flow rate : 1.0 mL/min.

Col. Temp. : Ambient

Detection : FL Ex 300 nm

Em 330 nm

Injection Vol. :

Sample : Tocopherols

Solutes : 1. DL- α - Tocopherol
2. D- β - Tocopherol
3. D- γ - Tocopherol
4. D- δ - Tocopherol