


Fat-soluble Vitamins

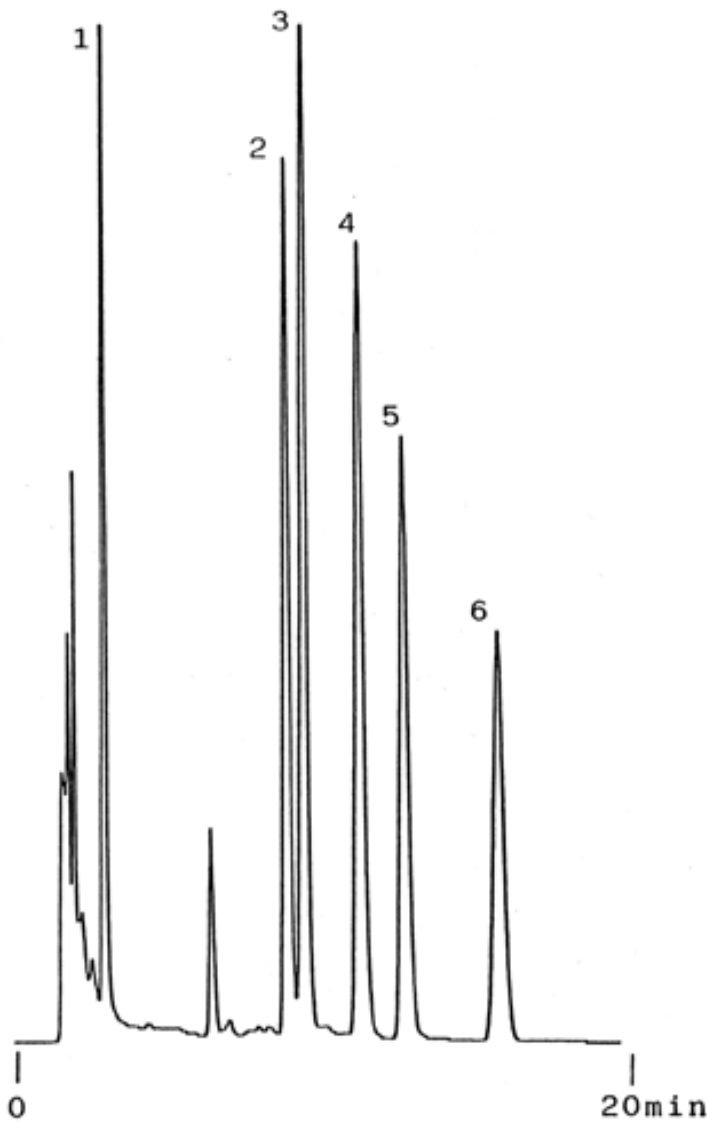

Column : Inertsil ODS-2
 5µm
 150 × 0.7 mm I.D.

Column Cat. No. : 5020-13073

Eluent : CH₃CN/CH₃OH = 75/25

Flow rate : 30 µL/min.

Col. Temp. :

Detection : UV 280 nm

Injection Vol. :

Sample : Fat-soluble Vitamins

Solutes : 1. Vitamin A
 2. Vitamin D2
 3. Vitamin D3
 4. Vitamin E
 5. Vitamin E acetate
 6. Vitamin K1