


Tricyclic antidepressant

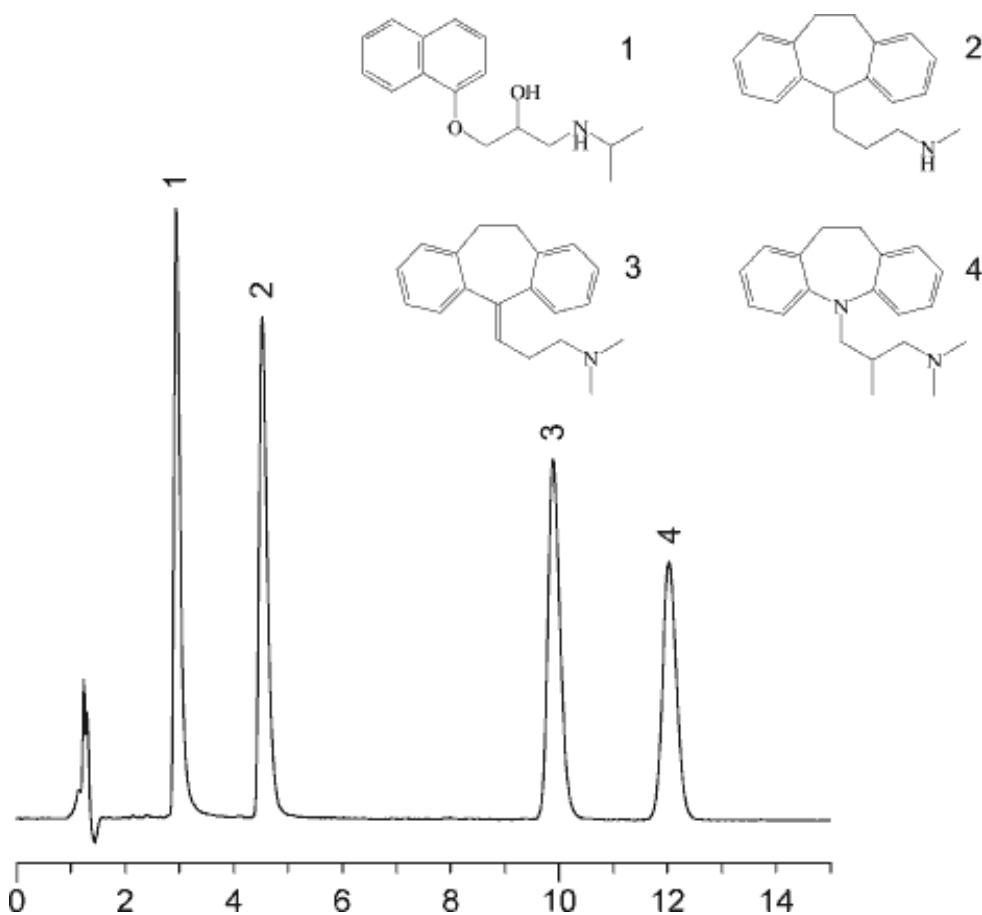

Column :Inertsil ODS-3
5µm
150 × 4.6mmI.D.

Column Cat. No. :
Eluent :20mM K₂HPO₄ (pH7.0) / Acetonitrile (30:70)
Flow rate :1.0 mL/min.
Col. Temp. :40°C
Detection :UV 254 nm
Injection Vol. :3µL
Sample :1. Propranolol (56µg/mL)

Solutes : 1. Propranolol (56µg/mL)
2. Nortriptyline (11µg/mL)
3. Amitriptyline (11µg/mL)
4. Trimipramine (11µg/mL)

2. Nortriptyline (11µg/mL)
3. Amitriptyline (11µg/mL)
4. Trimipramine (11µg/mL)