


Analysis of Fucoxanthin in Food

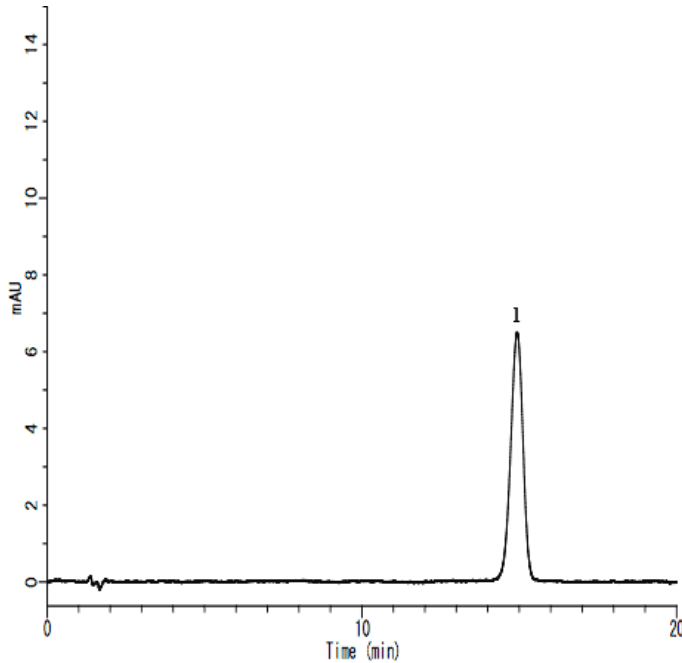

System : GL-7400 system
Column : Inertsil ODS-3 (4µm, 150 x 4.6 mm I.D.)
Column Cat. No. :
Eluent : A) CH₃CN
 : B) H₂O
 : A/B = 75/25, v/v
Flow rate : 1.0 mL/min.
Col. Temp. : 40 °C
Detection : PDA 445 nm (GL-7452A PDA Detector)
Injection Vol. : 20 µL
Sample :

Solutes : 1. Fucoxanthin (1mg/L)