


## Analysis of 8 kinds of Pesticides

Analyte: :

- 1 2,4-D
  - 2 Bentazon
  - 3 Triclopyr
  - 4 Bensulfuron methyl
  - 5 Flazasulfuron
  - 6 Mecoprop (MCP)
  - 7 Halosulfuron methyl
  - 8 Bensulide
- each 100 µg/L

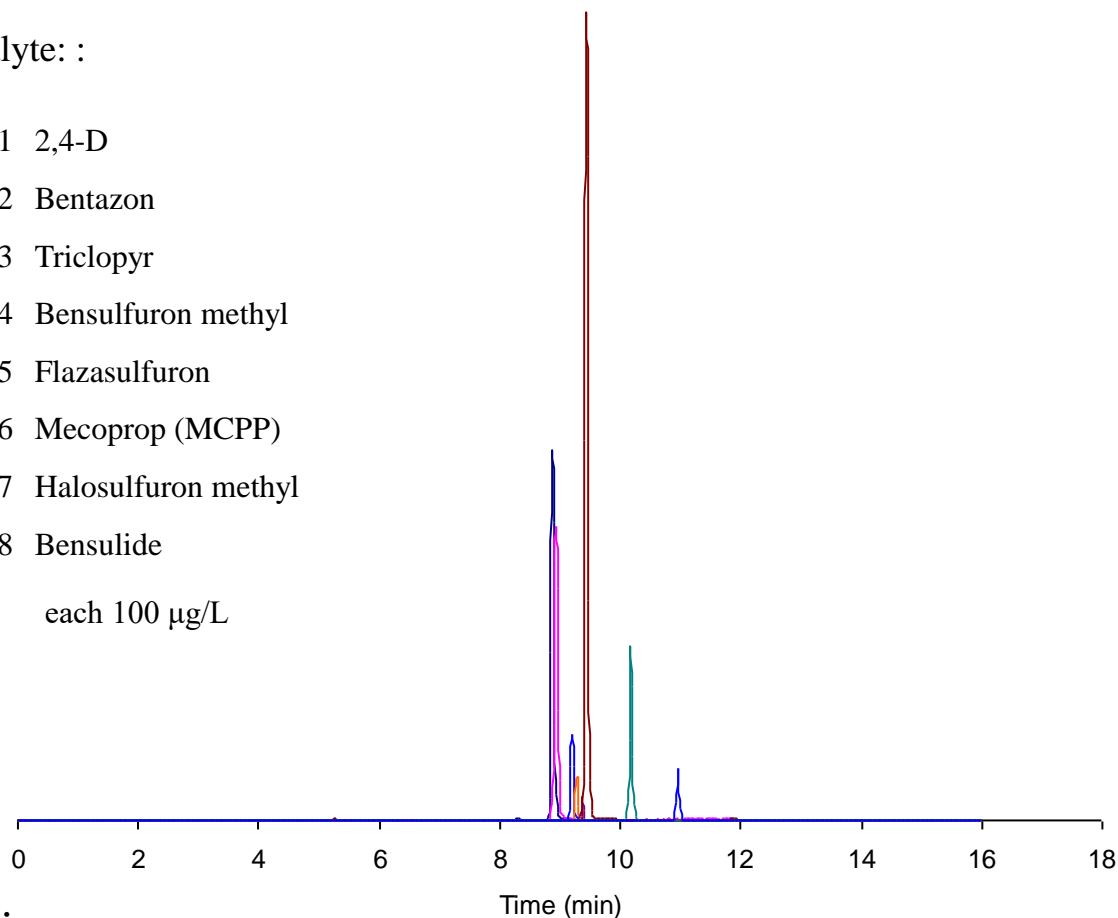

### Conditions

**System** : Symbiosis Pharma system  
**Column** : Inertsil ODS-SP (3 µm, 150 x 2.1 mm I.D)  
**Column Cat. No.** : 5020-02715  
**Eluent** : A) CH<sub>3</sub>CN  
           : B) 0.1% CH<sub>3</sub>COOH  
           : A/B = 0/100 - 10 min. - 100/0 - 5 min. hold, v/v  
**Flow Rate** : 0.2 mL/min  
**Col. Temp.** : 40 °C  
**Detection** : LC/MS/MS (ESI, Negative., MRM, Q1/Q3 conditions were in Supplemental Data)  
**Injection Vol.** : 10 µL  
**Analyte** : Pesticides

# Analysis of 8 kinds of Pesticides

## Supplemental Data

